


ENDOVASCULAR

brite tip[®]
S H E A T H

Interventional Sheath Introducer

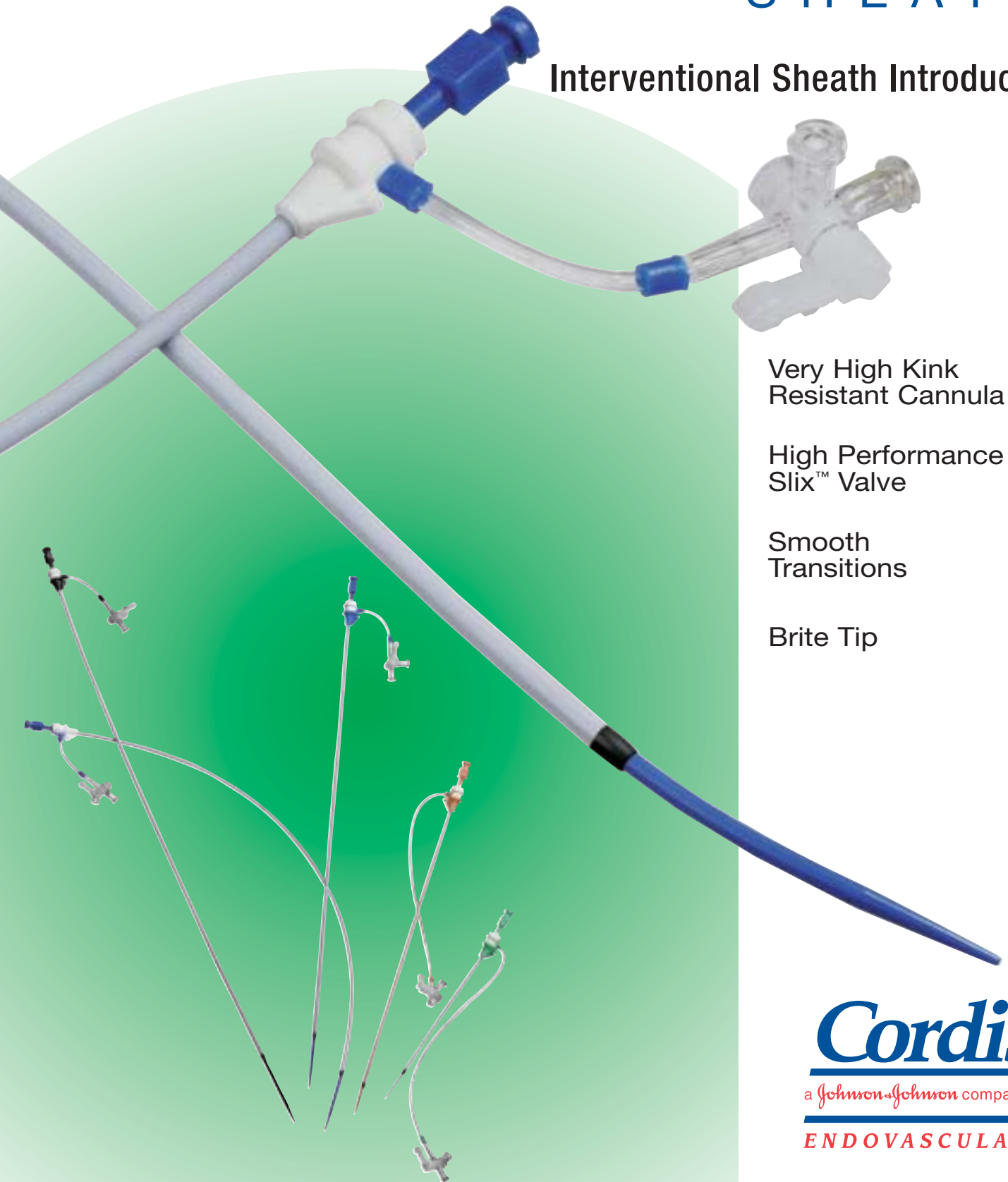

Very High Kink
Resistant Cannula

High Performance
Slix™ Valve

Smooth
Transitions


Brite Tip

Cordis[®]
a Johnson & Johnson company
ENDOVASCULAR

brite tip[®]

SHEATH INT

*designed to be
used for
interventional*

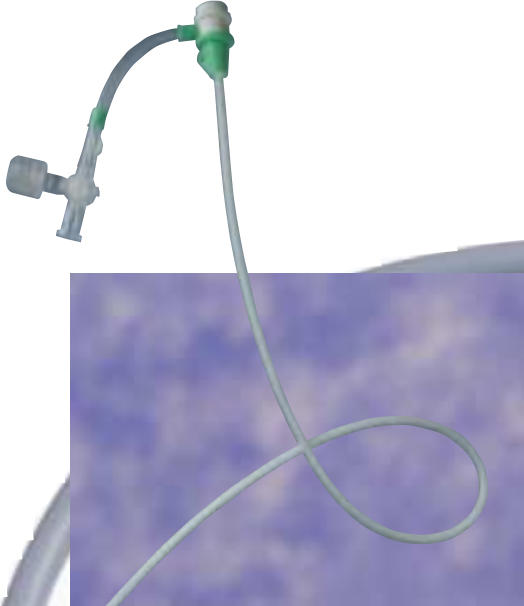

KINK RESISTANT CANNULA

Co-extruded (2 layers) cannula

- ▶ facilitates intravascular device passage, drug infusion and blood sampling

MDX coated cannula* and vessel dilator**

- ▶ enhances lubricity during insertion and provides anti-thrombogenic surface

*inside and outside
**outside

HIGH PERFORM



SMOOTH TRANSITIONS

Smooth cannula-to-dilator and dilator-to-guidewire transition

- ▶ provides effortless vessel entry and results in reduced groin complications

MULTIPLE

PRODUCER

*primarily
procedures*

APPLICATIONS

- Interventional procedures (PTA, stenting)
- Antegrade punctures
- Obese patients
- Patients with tortuous iliac arteries
- Trans-sheath contrast injection for verification of stent position (e.g. renal stenting)

ANCE SLIX™ VALVE

Six-cut spiral design (tailored to French size)

▶ minimises blood loss

Softer, more pliable silicone material

▶ is more lubricious

Minimal force required for inserting, withdrawing or turning a catheter

▶ provides superior manoeuvrability for smooth catheter placement

BRITE TIP™

3 mm atraumatic, radiopaque distal tip

▶ avoids positioning of intravascular device inside the sheath cannula








































LENGTHS PER FRENCH SIZE

Available French sizes : 4, 5, 6, 7, 8, 9, 10 and 11 F

Range of lengths : 5.5, 11, 23, 35, 45, 55, 90 cm

BRITE TIP INTERVENTIONAL SHEATH INTRODUCERS

Catalogue Numbers

French Size	Colour Code	Cannula usable length (cm)	Dilator Tip length* (mm)	Mini-Guidewire	Units/Package	
6		5.5	35	yes	5	401-605 M
7		5.5	35	yes	5	401-705 M
8		5.5	35	yes	5	401-805 M
<hr/>						
4		11	25	yes	5	401-411 M
5		11	25	yes	5	401-511 M
6		11	35	yes	5	401-611 M
7		11	35	yes	5	401-711 M
8		11	35	yes	5	401-811 M
9		11	35	yes	5	401-911 M
10		11	45	yes	5	401-011 M
11		11	45	yes	5	401-111 M
<hr/>						
4		23	25	no	5	401-423 M
5		23	25	no	5	401-523 M
6		23	35	no	5	401-623 M
7		23	35	no	5	401-723 M
8		23	35	no	5	401-823 M
9		23	35	no	5	401-923 M
10		23	45	no	5	401-023 M
11		23	45	no	5	401-123 M
<hr/>						
5		35	25	no	2	401-535 M
6		35	35	no	2	401-635 M
7		35	35	no	2	401-735 M
8		35	35	no	2	401-835 M
9		35	35	no	2	401-935 M
<hr/>						
5		45	25	no	2	401-545 M
6		45	35	no	2	401-645 M
7		45	35	no	2	401-745 M
8		45	35	no	2	401-845 M
9		45	35	no	2	401-945 M
<hr/>						
5		55	25	no	2	401-555 M
6		55	35	no	2	401-655 M
7		55	35	no	2	401-755 M
8		55	35	no	2	401-855 M
9		55	35	no	2	401-955 M
<hr/>						
6		90	35	no	2	401-690 M
7		90	35	no	2	401-790 M
8		90	35	no	2	401-890 M

Guidewire Acceptance: 0.035" * Length of the dilator protruding distally from the cannula.

Important Information: Prior to use, refer to the "Instructions for use" supplied with these devices for indications, contraindications, side effects, suggested procedure, warnings and precautions.
As part of the Cordis policy of continuous product development we reserve the right to change product specifications without prior notification.

Europe
Tel: +32 2 352 14 11
Fax: +32 2 352 15 90